

AD _____

Award Number: DAMD17-02-1-0386

TITLE: Genome Wide Examination of Allelic Loss in Lobular and
Ductal Breast Cancer

PRINCIPAL INVESTIGATOR: Peggy Porter, M.D.

CONTRACTING ORGANIZATION: Fred Hutchinson Cancer Research Center
Seattle, WA 98109-1024

REPORT DATE: July 2003

TYPE OF REPORT: Annual

PREPARED FOR: U.S. Army Medical Research and Materiel Command
Fort Detrick, Maryland 21702-5012

DISTRIBUTION STATEMENT: Approved for Public Release;
Distribution Unlimited

The views, opinions and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy or decision unless so designated by other documentation.

20040706 092

REPORT DOCUMENTATION PAGE

Form Approved
OMB No. 074-0188

Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503

1. AGENCY USE ONLY
(Leave blank)

2. REPORT DATE
July 2003

3. REPORT TYPE AND DATES COVERED
Annual (17 Jun 2002 - 16 Jun 2003)

4. TITLE AND SUBTITLE

Genome Wide Examination of Allelic Loss in Lobular and Ductal Breast Cancer

5. FUNDING NUMBERS

DAMD17-02-1-0386

6. AUTHOR(S)

Peggy Porter, M.D.

7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES)

Fred Hutchinson Cancer Research Center
Seattle, WA 98109-1024

E-Mail: PPORTER@FHCRC.ORG

8. PERFORMING ORGANIZATION
REPORT NUMBER

9. SPONSORING / MONITORING
AGENCY NAME(S) AND ADDRESS(ES)

U.S. Army Medical Research and Materiel Command
Fort Detrick, Maryland 21702-5012

10. SPONSORING / MONITORING
AGENCY REPORT NUMBER

11. SUPPLEMENTARY NOTES

12a. DISTRIBUTION / AVAILABILITY STATEMENT

Approved for Public Release; Distribution Unlimited

12b. DISTRIBUTION CODE

13. ABSTRACT (Maximum 200 Words)

The purpose of this study is to compare allelic loss throughout the genome in lobular and ductal type breast cancers using the commercially available HuSNP assay. We will evaluate the rate of allelic loss on individual chromosome arms in 60 each of lobular and ductal breast cancers obtained from a population-based study (NCI RO1 CA85913, P.I. Janet Daling). We have validated the prototype HuSNP platform for use in fixed tissue samples [Schubert, 2002 #1827], 2) enrolled and identified subjects for the study, 3) flow cytometry cell sorted breast tumors for the study, and 4) established the new array platform for collection of data from 10,000 SNPs in the Porter lab and the array facility at the FHCRC. Information about morphology-specific traits gained from this study of a large number of lobular cancers will lead to an increased understanding of the biology of distinct subsets of breast cancer. It will provide a basis for future studies designed to stratify patients into prognostic and treatment groups and/or inform the development of targeted therapies for specific tumor types.

14. SUBJECT TERMS

Invasive lobular carcinoma, invasive ductal carcinoma, loss of heterozygosity (LOH), post menopausal breast cancer, genomic analysis, genetic biomarkers

15. NUMBER OF PAGES
17

16. PRICE CODE

17. SECURITY CLASSIFICATION
OF REPORT
Unclassified

18. SECURITY CLASSIFICATION
OF THIS PAGE
Unclassified

19. SECURITY CLASSIFICATION
OF ABSTRACT
Unclassified

20. LIMITATION OF ABSTRACT
Unlimited

NSN 7540-01-280-5500

Standard Form 298 (Rev. 2-89)
Prescribed by ANSI Std. Z39-18
298-102

Table of Contents

Front Cover.....	1
SF 298.....	2
Table of Contents.....	3
Introduction.....	4
Body.....	4-5
Key Research Accomplishments.....	5
Reportable Outcomes.....	5
Conclusions.....	5
References.....	5
Appendices.....	6

Procedure for identification of eligible subjects, specimen transfer, and DNA extraction of buffy coats collected as part of the SHARE study.

Publication:

Schubert, E.L., Hsu, L., Cousens, L.A., Glogovac, J., Self, S., Reid, B.J., Rabinovitch, P.S. and Porter, P.L. Single Nucleotide Polymorphism Array Analysis of Flow-Sorted Epithelial Cells From Frozen *Versus* Fixed Tissues For Whole Genome Analysis of Allelic Loss in Breast Cancer. American Journal of Pathology, 160:73:79, 2002.

INTRODUCTION

Invasive lobular cancer (ILC) differs both morphologically and clinically from the more common ductal (IDC) type of breast cancer. ILC is more likely to be steroid hormone receptor positive, to have a low proliferation rate, to be diploid and to exhibit loss of e-cadherin protein expression. Loss of function of tumor suppressor genes or level of genomic instability may differ as well, but studies to date have examined relatively few lobular cases at only limited regions of the genome. In this study we are conducting a genome-wide assessment of allelic loss using a microarray (Affymetrix 10K SNP assay) composed of 10,000 single nucleotide polymorphisms (SNPs) to evaluate the rate of allelic loss on individual chromosome arms in 60 each of lobular and ductal breast cancers obtained from a population-based study (NCI R01 CA85913, P.I. Janet Daling). Data from the parent study will allow us to correlate allelic loss findings with risk factors such as, hormone replacement therapy (HRT) use, and other tumor characteristics such as estrogen and/or progesterone receptor status. A comparison of genome-wide allelic loss in ILC and IDC will provide new information about regions of genetic loss and genomic instability associated with particular phenotypic traits. Information about morphology-specific traits gained from studying a large number of lobular cancers will lead to an increased understanding of the biology of distinct subsets of breast cancer, and provide a basis for future studies that would define patient stratification into prognostic and treatment groups and/or inform the development of targeted therapies for specific tumor types.

BODY

The following describes the research accomplishments associated with each item in the approved Statement of Work.

Task 1: Selection of 60 lobular and 60 ductal cases, (months 1-14)

We have compiled data on women in the parent study as they are enrolled. We use an algorithm that selects subjects based on number of flow cytometrically sorted tumor cells available, estrogen receptor status of tumor, stage of disease and age (see appendix for protocol). To date, 44 ILC and 45 IDC tumors have been identified. DNA has been extracted from the tumors and from the blood lymphocytes for the subject.

Task 2: Prepare DNA for HuSNP™ analysis (months 6-22)

DNA is stored at -80°C until all samples are collected. Whole-genome amplification (PEP) will be performed in batches when all 120 samples are available.

Task 3: HuSNP™ analysis months (9-22)

The performance of the HuSNP assays and data collection has been delayed until the new Affymetrix array with 10,000 SNPs is brought on line. Although the higher density array was scheduled for commercialization in early 2002, production was delayed and the company is now establishing the array hybridization protocols and software in institutions around the U.S. Dr. Loo and technicians in the DNA array facility at the FHCRC will be trained on the new array

assays November 18th through 20th and we anticipate assay availability as of December 2003. One hundred 10K arrays have been purchased and are ready for assays when the technology is fully tested in the facility.

Task 4. Data analysis (months 9-24)

Data analysis has not begun.

Task 5. Publication of results (months 22-24)

No publications have resulted from the project at this time.

KEY RESEARCH ACCOMPLISHMENTS

- validation of the prototype HuSNP platform for use in fixed tissue samples¹
- validation in progress for the new higher density 10,000 SNP array
- enrollment and identification of subjects for study using 10,000 SNP arrays
- flow cytometry cell sorting of tumors for study
- establishment of new array platform for collection of data from 10,000 SNPs in the Porter lab and the array facility at the FHCRC

REPORTABLE OUTCOMES

Schubert EL, Malone K, Daling JD, Cousens LG, Porter PL. Whole genome LOH analysis of lobular and ductal breast cancers by Hu-SNPTM array. Poster 24th Annual San Antonio Breast Cancer Symposium, 2001.

Schubert EL, Hsu L, Cousens LA, Glogovac J, Self S, Reid BJ, Rabinovitch P, Porter PL. Single nucleotide polymorphism array analysis of flow-sorted epithelial cells from frozen *versus* fixed tissues for whole genome analysis of allelic loss in breast cancer. *Am J Pathol*, 160(1): 73-79, 2002.

CONCLUSIONS

With the establishment of the array platform at the FHCRC, the identification of subjects and the processing of DNA from selected tumors, we are close to the data collection phase of the study. We will generate new data concerning the molecular changes associated with lobular and ductal breast cancer that will lead to an increased understanding of the biologic differences between these subsets. Most optimistically, our findings could inform the development of targeted detection and therapy strategies for specific tumor types.

REFERENCES

1. Schubert EL, Hsu L, Cousens LA, et al. Single Nucleotide Polymorphism Array Analysis of Flow-Sorted Epithelial Cells from Frozen Versus Fixed Tissues for Whole Genome Analysis of Allelic Loss in Breast Cancer. *Am J Pathol*. 2002;160(1):73-79.

APPENDICES

Procedure for identification of eligible subjects, specimen transfer, and DNA extraction of buffy coats collected as part of the SHARE study.

Schubert EL, Hsu L, Cousens LA, et al. Single Nucleotide Polymorphism Array Analysis of Flow-Sorted Epithelial Cells from Frozen Versus Fixed Tissues for Whole Genome Analysis of Allelic Loss in Breast Cancer. *Am J Pathol.* 2002;160(1):73-79.

SHARE / LOBULAR HuSNP PROJECT

Procedure for IDENTIFICATION OF ELIGIBLE SUBJECTS, SPECIMEN TRANSFER, AND DNA EXTRACTION OF BUFFY COATS COLLECTED AS PART OF THE SHARE STUDY

All tissue requests and blood collections for the SHARE study originate with CERC Studies staff. The Porter Lab will receive the Formalin-fixed, paraffin-embedded tumor tissue in the lab and test DNA extracted from flow-sorted tumor cells against a normal (DNA extracted from the blood lymphocytes by the Specimen Processing Lab) for each of 120 eligible subjects (60 invasive ductal carcinomas (8500/3), and 60 invasive lobular carcinomas (8520/3) to be tested by HuSNP array. The selection criteria for inclusion into the HuSNP project are: enrollment in SHARE study, unilateral breast cancer, diagnosis date between January 1 2000 and June 30 2001, Stage II or lower, ER positive, minimum 50K flow-sorted tumor cells available, and blood sample collected and stored by CERC.

1. CASE SELECTION / SUBJECT IDENTIFICATION

a. Identification of potential subjects by the Porter Lab

Tumor tissue is received by the Porter Lab. After histological review of each tumor, IHC and flow cytometry, the Porter Lab identifies ER positive, pure lobular and pure ductal invasive breast tumors from the SHARE cohort for which a minimum of 50,000 sorted cells have been stored. Approximately once every 3 months, the Porter Lab prepares a list of the tumors meeting these criteria, identified by Study ID along with the correlating histological dx (per Porter Lab pathologists) for each case. The list is sent via courier, on a disc, to CERC Programmer.

b. Selection of subjects by CERC

The potential subjects identified by the Porter lab are screened by CERC programmer for additional eligibility parameters (dx date, stage, laterality, and current blood sample availability). He will also consider the histology code sent by the Porter Lab, in the interest of acquiring 60 suitable cases of each diagnosis (invasive ductal and invasive lobular). Patients eliminated by any criterion, or any who have refused blood donation, are noted as permanently ineligible. Subjects from whom blood collection has not yet occurred are noted as temporarily ineligible. Any of these for whom blood is later collected are included with the next group of specimens, after CERC Study Manager notifies CERC Programmer, who changes their status to "eligible" and flags them for inclusion in the next group.

c. Preparation of labels and tracking lists

- CERC Programmer

- (1) sends CERC Study Manager an email notifying her of any "temporarily ineligible" subjects for blood collection follow-up. He also sends this information to Porter Lab.
 - (2) generates a list of the Study IDs and Buffy Coat IDs of the eligible participants requested, which he downloads onto a disc and sends via courier to Porter Lab.
 - (3) emails a "Buffy Coat Pack List" to CERC Study Manager listing the Buffy Coat ID's of the specimens to be pulled and sent to the Specimen Processing Lab.
 - (4) prints 2 identical tube labels for each blood specimen, on which is printed the 5 digit "Buffy Coat ID" number (with which the buffy coat is already also labeled) and a CERC "DNA Aliquot ID" number, which has been previously set aside for this purpose.
- CERC Data Coordinator
 - (1) prints several sheets of duplicate sequential "DNA Aliquot ID" labels. The ID numbers are prefaced by the letters "SHA". Half of these sheets of labels are printed in bold red type with the suffix "QIA" following the ID number on the label. The other labels are printed in standard black type with the suffix "PHE" following the number on the label.
 - (2) prints the "SHARE DNA Tracking Log Forms", preparing one for each Buffy Coat ID, pre-printing the DNA Aliquot ID to appear on the Porter Lab's 5 µg aliquot on the form in the "DNA Aliquot ID for Porter Lab" field, as well as the Buffy Coat ID.
 - Cancer Biology Division Programmer downloads the list of eligible subjects (received from CERC Programmer) into the Porter Lab SHARE database.
 - Porter Lab Data Coordinator
 - (1) picks up the disc (if any) from the courier and delivers it to Stephanie Stafford.
 - (2) prepares 5 tube labels with the Porter Lab's "Lab Number", scheduled extraction and aliquotting date, and the correlating Buffy Coat ID, for each eligible subject listed.
 - CERC Study Manager
 - (1) attaches the two special labels which match each Buffy Coat ID to each correlating "SHARE DNA Tracking Log Form" with a paperclip.
 - (2) pulls the Buffy Coats listed on the "Buffy Coat Pack List" from the freezer and packs them carefully in dry ice in a cooler as described below.
 - (3) sends the Forms and Aliquot ID labels, specimens and Pack List to the Specimen Processing Lab as described below

2. SPECIMEN TRANSFER

Study Manager contacts Porter Lab via phone or email after receiving the "Buffy Coat Pack List" to let the lab know exactly when to send the blood specimens to the Specimen Processing Lab. Porter Lab then contacts the Specimen Processing Lab to schedule the

extractions at a time when Porter Lab staff can pick up Porter Lab aliquot immediately after extraction. Study Manager then pulls the appropriate blood specimens from the freezer and packs them in a cooler in dry ice. She includes the forms and labels prepared by CERC programmer in the package and labels everything clearly with "VIA COURIER, To: The Specimen Processing Lab, DE-120, Day Campus, x 4645, Specimens Enclosed, Keep Frozen, Contact: SHARE, Porter Lab, x5470". The package is then left with the receptionist at the Met Park 2nd floor Reception Desk for the next courier run. Alternatively, the specimens, forms, and labels are delivered to the Specimen Processing Lab by CERC Staff.

3. EXTRACTION / ELUSION / LAB TRACKING

The Specimen Processing Lab will:

- Perform the extractions at the scheduled time, notifying the Porter Lab of any schedule changes by calling Porter Lab at 667-5470 or 667-6187.
- Extract one quarter of each buffy coat via the Qiagen method, and three quarters by hand using the phenol chloroform method.
- Indicate the date of the extraction, the name of the Specimen Processing Lab Tech performing the extraction, the Elution used, and total volume and total μg extracted, and the final concentration on the "SHARE DNA Tracking Log Form" for each buffy coat, for each extraction performed.
- Perform the extractions, and aliquot all resulting material from the Qiagen method into 5 μg aliquots and all resulting material from the Phenol chloroform method into Master and Children aliquots according to that protocol.
- Phone Porter Lab at 667-5470 or 667-6187 to notify them that the extractions are complete
- Ensure that the first aliquot of DNA purified by the Qiagen method from each Buffy Coat contains a total of 5 μg of DNA where possible, at a concentration of 50 $\text{ng}/\mu\text{l}$, designating this aliquot for the Porter Lab. The Porter Lab aliquot remains UNFROZEN.
- Label this first aliquot of DNA purified by the Qiagen method from each Buffy Coat for the Porter Lab, using the provided label (attached to the correlating "SHARE DNA Tracking Log Form" with a paperclip), which has both the originating Buffy Coat ID number and a regular "DNA Aliquot ID" printed on it. One such label goes onto the microcentrifuge tube, and the other matching label with the same DNA Aliquot ID goes onto the "SHARE DNA Tracking Log Form".
- Label all other aliquots with a standard "DNA Aliquot ID" label (provided), which begins with the letters "SHA". One such label is stuck onto the tube, and the other

matching label with the same DNA Aliquot ID is stuck onto the "SHARE DNA Tracking Log Form" on the line for that aliquot.

- Complete additional fields on the "SHARE DNA Tracking Log Form" ("DNA Amt", "DNA Conc", "Aliquot Vol", and "Box Name or number") for each aliquot.
- Temporarily store the Porter Lab's aliquot in a separate box marked "SHARE Qiagen Porter" at room temperature until Porter Lab staff arrives to perform further aliquotting. DO NOT FREEZE THE PORTER LAB'S ALIQUOT.
- Hand off the Porter Lab's aliquot to Porter Lab staff for further aliquotting.
- Store remaining Qiagen method aliquots in the freezer, in a box labeled "SHARE Qiagen W-D", and all aliquots resulting from the phenol chloroform method extraction in a freezer box labeled "SHARE Phenol W-D".
- Notify CERC Study Manager (667-5044) when CERC's tubes are ready for pick-up and transfer to their own freezers for long-term storage.

4. RECEIPT / STORAGE

- a) Porter Lab Research Tech will be present when the extractions take place so that the Porter Lab's 5- μ g aliquot for each participant can be picked up at room temperature and re-aliquotted by them into 5 separate microcentrifuge tubes, which they then label with the DNA Aliquot ID, concentration (ng/ μ g), current date, and the Porter Lab's "Lab number".
- b) When aliquotting of all extracted material is finished, the completed "SHARE DNA Tracking Log Forms" and the Porter Lab's aliquots are taken back to the Porter Lab for "check in" and storage. Discrepancies between the imported list of Buffy Coat ID's and those appearing on the Forms are noted and resolved after investigation by Porter Lab and CERC Study Manager, and records are updated accordingly if necessary. The "SHARE DNA Tracking Log Forms" are photocopied, and the copies are filed in the Porter Lab by date of extraction. The original Forms are sent to CERC Study Manager for data entry and storage.
- c) The two boxes containing CERC's aliquots, which are temporarily stored in the Specimen Processing Lab's freezers, are picked up by CERC staff and transported by hand back to their own freezers. They are then put into long term storage there while testing methods are being developed for further study, and for archival purposes in case additional DNA is needed by the Porter Lab for testing.

Technical Advance

Single Nucleotide Polymorphism Array Analysis of Flow-Sorted Epithelial Cells from Frozen Versus Fixed Tissues for Whole Genome Analysis of Allelic Loss in Breast Cancer

Elizabeth L. Schubert,* Li Hsu,[†]
Laura A. Cousens,* Jeri Glogovac,* Steve Self,[†]
Brian J. Reid,*[‡] Peter S. Rabinovitch,*[§] and
Peggy L. Porter*[§]

From the Divisions of Human Biology and Public Health Sciences, Program in Cancer Biology, and the Division of Public Health Sciences,[†] Program in Biostatistics, Fred Hutchinson Cancer Research Center, Seattle, and the Departments of Gastroenterology[‡] and Pathology,[§] University of Washington, Seattle, Washington*

Analysis of allelic loss in archival tumor specimens is constrained by quality and quantity of tissue and by technical limitations on the number of chromosomal sites that can be efficiently evaluated in conventional analyses using polymorphic microsatellite markers. Newly developed array-based assays have the potential to yield genome-wide data from small amounts of tissue but have not been validated for use with routinely processed specimens. We used the Affymetrix HuSNP assay, composed of 1494 single nucleotide polymorphism sites, to compare allelic loss results obtained from both formalin-fixed and frozen breast tissue samples. Tumor cells were separated from normal epithelia and nonepithelial cells by dissection and bivariate cytokeratin/DNA flow sorting; normal breast cells from the same patient served as constitutive normal. Allele results from the HuSNP array averaged 96% reproducibility between duplicates and were concordant between the fixed and frozen normal samples. We also analyzed DNA from the same samples after whole-genome amplification (primer extension preamplification). Although overall signal intensities were lower, the genotype data from the primer extension preamplification material was concordant with genomic DNA data from the same samples. Results from genomic normal tissue DNA averaged informative single nucleotide polymorphism at 379 (25%) loci genome-wide. Although data points

were clustered and some segments of chromosomes were not informative, our data indicated that the Affymetrix HuSNP assay could provide an efficient and valid genome-wide analysis of allelic imbalance in routinely processed and whole genome-amplified pathology specimens. (*Am J Pathol* 2002, 160:73-79)

Loss of heterozygosity (LOH), or allelic loss, is one of the most frequent genetic abnormalities in breast cancer. It may serve as a marker of generalized genomic instability, and when frequently observed in a region, it is considered indirect evidence for the presence of a tumor suppressor gene within that region of loss. In sporadic breast cancer, allelic loss at multiple chromosomal locations has been identified in a range of invasive and preinvasive breast cancers as well as benign and normal breast epithelium adjacent to tumor.¹⁻³ However, a complete evaluation of LOH in breast cancer has been hampered by the limited number of polymorphic markers available for study; the heterogeneity of breast tissue (mixed non-tumor and tumor cells); the lack of sufficient numbers of fresh or frozen samples with associated demographic or clinical data, and the small amount of tissue available from currently diagnosed breast cancers. To address these limitations, we used flow cytometry to select and purify tumor cells from routinely processed tissue blocks, whole genome amplification to increase the amount of DNA available for study, and a microarray assay to assess all chromosomes efficiently and simultaneously.

The newly developed Affymetrix HuSNP array, which contains 1494 single nucleotide polymorphism (SNP)

Supported by NIH/NCI grants RO1 CA78855 (to P. S. R.) and RO1 CA71735 (to P. L. P.).

Accepted for publication September 25, 2001.

Address reprint requests to Dr. Peggy L. Porter, Department of Pathology, Box 357705, University of Washington, 1959 NW Pacific St., Seattle, WA 98195-7705. E-mail: pporter@fhcrc.org.

sites genome-wide and requires only 135 ng of genomic DNA (gDNA) per assay, is a potential platform for evaluating genome-wide genetic analysis of breast tissue. The usefulness of a prototype SNP array and the current HuSNP array for analysis of allelic loss in fresh lung tumors removed at autopsy and fresh biopsies from esophageal cancers, respectively, has been previously described.^{4,5} However, the analysis of formalin-fixed, paraffin-embedded pathology specimens by the commercially available HuSNP assay has not been reported.

Here we discuss the use of the HuSNP to examine allelic imbalance in both frozen and fixed pathology specimens and compare results between the two preservation methods. To purify populations of cells from the tissue for analysis we used bivariate flow cytometry, which allowed us to sort tumor cells for analysis based on positive cytokeratin staining and gDNA content.⁶ In addition to gDNA, we also examined the use of a polymerase chain reaction (PCR)-based whole-genome amplification method, primer-extension-preamplification (PEP) that increases the amount of template available for analysis ~30-fold^{7,8} and compared allelic loss results from the PEP product to results with gDNA. HuSNP allelic loss results were also compared to results from conventional polymorphic microsatellite markers (short tandem repeats or STRs) on chromosomes 11 and 17.

Materials and Methods

Tissue Samples

Tumor and normal tissue from two breast cancer patients were obtained from the University of Washington tissue bank with patient consent and in compliance with the Institutional Review Board. Samples taken at the time of surgery were divided into two portions and each portion was processed routinely either by freezing in OCT media or formalin fixation followed by paraffin embedding. No gross difference was apparent between the portions selected for either preservation method. The presence of tumor in each block was confirmed microscopically. A formalin-fixed tissue block from each of the cases was tested for estrogen receptor, progesterone receptor, c-erbB2 oncogene protein, and p53 tumor suppressor gene protein by immunohistochemistry as previously described.⁹

Flow Cytometry

Flow cytometry was performed on the frozen and fixed samples to purify tumor cells from normal epithelia and nonepithelial cells. Hematoxylin and eosin (H&E)-stained slides from both frozen and paraffin-embedded tumor sections were examined to confirm that the samples contained tumor epithelium. Similarly, H&E slides taken from the normal block confirmed that the sample contained no tumor.

From each frozen breast tissue sample, 20 to 50 50- μ m sections were cut and placed into phosphate-

buffered saline containing 1% bovine serum albumin (PBA).¹⁰

The samples were mechanically disaggregated and washed in PBA. The resulting cell suspensions were fixed in 0.5% electron microscopy grade formaldehyde and permeabilized in 0.1% triton/PBA before staining. From formalin-fixed tissue blocks, flow cytometry preparation was performed as described.⁶ Briefly, 1 to 20 60- μ m sections were cut from normal and tumor tissue blocks; regions of tumor in each section were dissected from surrounding tissue with a scalpel blade. All sections were deparaffinized, rehydrated, and digested in collagenase before a brief pepsin digestion.

Cell suspensions from both the frozen and fixed samples were stained with 4,6-diamidino-2-phenylindole and R-phyco-erythrin labeled AE1/AE3 (Roche, Indianapolis, IN), which recognizes a wide variety of acidic and basic cytokeratins. A parallel sample of cells was stained with R-PE-labeled isotype-matched mouse Ig (R-PE labeled IgG1; DAKO, Carpinteria, CA) and used as a negative control. Before sorting, all samples were forced through a 25-gauge needle 10 times to ensure a single cell suspension.

Cytokeratin-positive tumor cells were sorted by bivariate analysis with 488 nm and UV excitation on a Becton Dickinson (Mountain View, CA) FACS Vantage. R-PE, cytokeratin-positive populations were sorted based on their 4,6-diamidino-2-phenylindole-fluorescent DNA content, expressed as DNA index (DI = mean aneuploid G₁ fluorescence/mean diploid G₁ fluorescence). Cells from the normal blocks were processed and stained similarly to the tumor samples. The DNA from all cells in the normal blocks was used as the constitutive normal for comparison with the tumor cell DNA.

Preparation of DNA Samples

DNA was extracted from frozen cells using the Puregene DNA isolation kit (Gentra Systems, Minneapolis, MN), following the manufacturer's suggestions with the addition of 1 μ l of 20 mg/ml of Proteinase K to the cell lysis buffer, followed by incubation at 50°C for 1 to 16 hours. DNA was extracted from fixed cells using a simple Proteinase K digestion method previously described.¹¹ Extracted DNA samples were quantified using the Picogreen dsDNA Quantitation Kit (Molecular Probes, Eugene, OR) on the Cytofluor II Fluorescence Multiwell Plate Reader (PerSeptive Biosystem Inc., Framingham, MA).

Whole genome amplification using the primer extension protocol (PEP) was performed as described.⁸ For each gDNA sample, six individual PEP reactions, each using 7 ng of gDNA as template, were performed and the PEP material pooled.¹² PEP material was used directly in the array protocol without purification or alteration of concentration.

STR Protocol

Twenty-four polymorphic repeat loci (STRs) on chromosomes 11 and 17 were amplified using fluorescent prim-

ers with PEP template. Chromosomes 11 and 17 were selected for allelic loss comparison between arrays and conventional repeat markers because both contain sites that are frequently lost in breast cancer.^{2,3} Markers were selected from those commercially available from Research Genetics (www.resgen.com) to obtain a survey of sites that corresponded as closely as possible to the HuSNP sites along the chromosomes. The physical locations of the markers in Mb are listed in Figure 1, B and C, as given by National Center for Biotechnology Information in July 2001 (www.ncbi.nlm.nih.gov). Primers for chromosome 11 were (11ptel) D11S1397, D11S2368, D11S2001, D11S1918, D11S1395, (cen), D11S4076, D11S1394, D11S4151, D11S2360, (11qtel). Those for chromosome 17 were (17ptel), D17S919, D17S1298, D17S1537, TP53, D17S786, D17S1541, D17S974, D17S975, (cen), D17S1293, D17S1158, D17S1294, D17S1185, D17S1305, D17S1290, D17S1288, (17qtel). PCR reactions were performed using standard protocols with PEP material as a template. PCR reaction products were multiplexed and then purified using Microcon-100 columns, after which the DNA was resuspended in sterile water. Reactions were run on an ABI 377 and analyzed using ABI Prism Gene Scan software.

To evaluate allelic loss for each marker^{12,13} the peak height of the first allele was divided by the peak height of the second allele to obtain the allelic ratio (AR). Samples were deemed informative at a locus if the AR for the normal tissue sample was sufficiently close to 1 (defined operationally as $0.7 < AR < 1/0.7$). For informative loci, an index Q was computed as the AR of the tumor tissue sample divided by the AR of the normal tissue sample.¹⁴ A locus was scored as having LOH if the Q value was sufficiently far from 1 (defined operationally as either $Q < 0.3$ or $Q > 1/0.3$). A locus was scored as retaining heterozygosity if the Q value was sufficiently close to 1 (defined operationally as $0.7 < Q < 1/0.7$).

HuSNP Protocol

The Affymetrix HuSNP protocol was performed according to manufacturer's instructions and as described.⁴ Each individual gDNA sample from both cases was analyzed twice in completely separate reactions, to yield data from a total of 18 HuSNP arrays for the two cases (five samples from case 1 in duplicate plus four samples from case 2 in duplicate). Similarly, the PEP material from each sample was analyzed by HuSNP in duplicate (18 HuSNP arrays). Data analysis using the Affymetrix Genechip software resulted in genotype calls that were used in the statistical analysis. The genetic map used in the analysis came from Affymetrix, release date June 2001.

Statistical Methods

To quantify the reproducibility of the HuSNP chips, the reliability measure was calculated. The reproducibility for making a consistent genotype call was defined as the number of SNPs with the same genotype calls from both replicates divided by the total number of SNPs for which

both replicates yielded signal calls. The reproducibility for making no-signal calls was also calculated, and was defined as the number of SNPs for which both replicates yielded no-signal calls divided by the total number of SNPs for which at least one replicate yielded a no-signal call.

Similarly, concordance of genotype and no-signal calls were measured between frozen and fixed tissue samples as well as between gDNA and PEP samples. Because each sample was analyzed in duplicate, there were a total of four possible comparisons between each set of fixed and frozen samples. The concordance measure was calculated by the ratio of the average over the four comparisons of the number of SNPs with same genotype calls from both samples and the average over the four comparisons of the number of SNPs for which both samples yielded signal calls. The concordance measure for no-signal calls was calculated similarly.

The informativity and allelic loss of the SNPs was examined for both cases. We defined a SNP as informative when one normal tissue replicate of the SNP was heterozygous (AB) and the other replicate was either heterozygous or had no signal. We defined a SNP site as having allelic loss when that SNP was informative in the normal tissue, one tumor tissue replicate of the SNP was hemizygous or homozygous (AA or BB), and the other replicate was either hemizygous, homozygous, or no signal.

All statistical analyses were performed using SPLUS statistical software (S-PLUS Reference Manual, version 3.2; Statistical Sciences I, Seattle, Washington).

Results

Subject and Tumor Characteristics

The results of the pathology review and immunohistochemical assays from the two cases used in this study are shown in Table 1. The patients' ages were similar at their respective times of diagnosis. Flow cytometric analysis revealed multiple aneuploid cell populations in the tumor from case 1. One cell population from the fixed tumor and one from the frozen tumor had very similar DIs (1.49 and 1.43, respectively), whereas an additional cell population seen only in the fixed portion of the tumor had a distinct DI of 1.82. Case 2 had a single aneuploid tumor

Table 1. Patient and Tumor Characteristics of the Two Breast Cancer Cases

	Case 1	Case 2
Histologic type	Lobular	Ductal
Age at diagnosis	48	46
AJCC stage	IIIA	IIA
ER protein status	Negative	Positive
PR protein status	Negative	Positive
c-erb-2 protein status	Negative	Positive
P53 protein status	Positive	Negative
DNA index	Multiple aneuploid	Aneuploid
DNA index—frozen tumor	1.43	1.79
DNA index—fixed tumor	1.49 and 1.82	1.76

ER, estrogen receptor; PR, progesterone receptor.

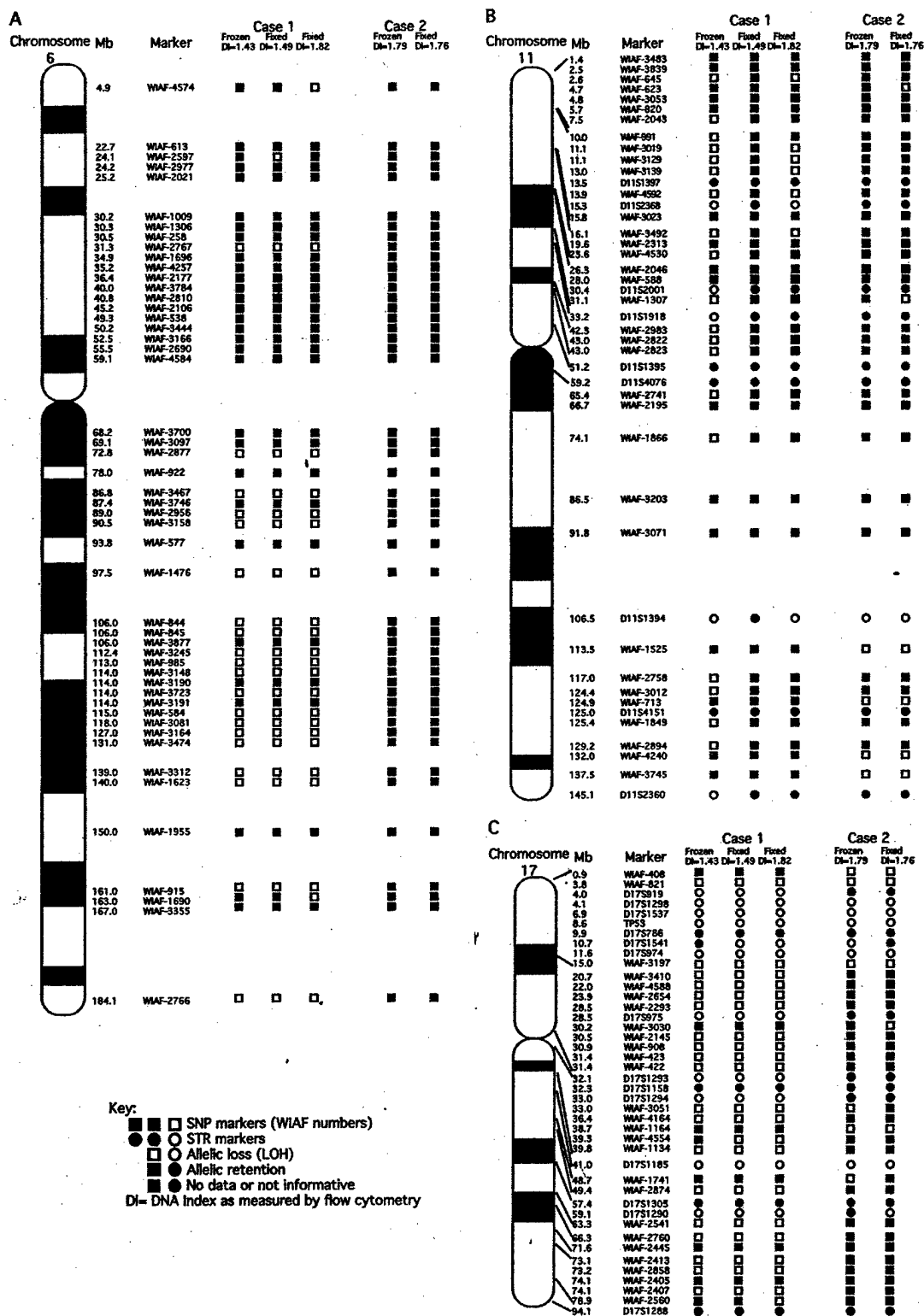


Figure 1. Graphical representation of HuSNP and STR allelic loss data for chromosomes 6 (A), 11 (B), and 17 (C) from both patients. SNP data are depicted as squares and STR data are indicated by circles. Shown are the data from all five cell populations isolated by bivariate flow cytometry in the two patients' tumors. SNP loci that were uninformative in both the frozen and fixed material from a patient were removed from this diagram. White squares and circles indicate allelic loss, black squares and circles indicate retention, and gray squares and circles indicate lack of informativity or data (in the STR markers) or no signal calls (in the HuSNP). Markers are arranged on the diagram to represent their physical distribution on the chromosome.

Table 2. Reproducibility of HuSNP Calls for Each Sample Type*

	Frozen normal		Frozen tumor		Fixed normal		Fixed tumor	
	Definite calls [†]	No signal calls	Definite calls	No signal calls	Definite calls	No signal calls	Definite calls	No signal calls
gDNA	98%	72%	96%	65%	98%	75%	95%	75%
PEP	97%	75%	95%	67%	95%	70%	93%	79%

*The average percent of definite (AA, AB, and BB) and no signal HuSNP genotype calls that were identical between independent, duplicate analyses of each sample type.

[†]AA, AB, or BB calls.

cell population distinguishable in both the fixed (DI = 1.76) and frozen (DI = 1.79) tissue samples. All three tumor cell populations from case 1 and both from case 2 were tested independently and included in the subsequent array analysis using both gDNA and PEP material from these cases.

HuSNP Analysis

For each SNP site on the chip, genotype results from the Affymetrix Genechip software were reported as definite calls (AA, AB, BB), no signal, or an intermediate call (AB.A or AB.B). The Genechip software does not score allele copy number but instead always indicates two alleles (AA and BB). The HuSNP chip contains 1494 individual SNP sites, however our experience was similar to that of a previous report,⁴ in that more than 100 of the 1494 sites on the chip consistently failed, yielding most of the no-signal calls. Intermediate calls were rare, seen in ~1% of sites in each assay. The reproducibility statistics for definite calls and no-signal calls between duplicate assays using the same gDNA or PEP sample are shown in Table 2. Table 3, A and B, show the concordance of definite and no-signal results between the fixed and frozen gDNA and PEP samples from each case (Table 3A) as well as the concordance between the gDNA and PEP results (Table 3B). A graphical representation of concor-

dance between fixed and frozen gDNA samples on chromosomes 6, 11, and 17 are presented in Figure 1.

The allelic loss results for each chromosome and gDNA sample source are shown in Table 4. PEP results were similar to the gDNA, as indicated by the genotype concordances shown in Table 3B. Informativity varied slightly between the fixed and frozen samples and resulted in some differences in allelic loss results between the paired samples as shown in Table 4. This variance was primarily because of no-signal calls at a particular SNP site in one sample type or the other and not to actual differences in calls between the fixed and frozen samples. The exception to this is in case 1 that contained more than one aneuploid population of tumor cells (Table 1). Allelic losses were reproducible in the duplicate analyses and generally concordant with adjacent sites in large regions along chromosomes, as is shown visually for chromosomes 6, 11, and 17 in Figure 1.

Comparison of HuSNP and STR Analyses

Of the 24 STR markers analyzed on chromosomes 11 and 17, there were a total of 69 informative sites between all five tumor populations identified in the two patients (see Figure 1, B and C, for details). Of these 69 sites, 60 showed correlation with data from adjacent HuSNP markers. However, at the nine STR sites that do not correlate

Table 3A. The Concordance[†] of Definite and No Signal Calls between the Duplicate Analyses of Frozen and Fixed Samples from Each of Two Breast Cancer Cases

	Case 1				Case 2			
	Normal tissue		Tumor tissue		Normal tissue		Tumor tissue	
	Definite calls*	No signal calls	Definite calls	No signal calls	Definite calls	No signal calls	Definite calls	No signal calls
gDNA	1162 (97%) [‡]	185 (64%)	1071 (92%)	199 (61%)	1150 (93%)	182 (71%)	1019 (95%)	189 (45%)
PEP	878 (92%)	543 (36%)	766 (87%)	288 (47%)	1131 (93%)	185 (68%)	919 (95%)	196 (37%)

Table 3B. Concordance[†] of gDNA and PEP DNA Data for Each Sample Type in Both Cases

	Frozen normal		Frozen tumor		Fixed normal		Fixed tumor	
	Definite calls*	No signal calls	Definite calls	No signal calls	Definite calls	No signal calls	Definite calls	No signal calls
Case 1	1004 (95%)	211 (49%)	1006 (92%)	219 (54%)	909 (91%)	188 (38%)	1200 (98%)	207 (78%)
Case 2	1183 (98%)	162 (57%)	1166 (97%)	167 (58%)	1166 (97%)	167 (58%)	825 (93%)	263 (43%)

*AA, AB, or BB calls.

[†]The concordance of genotype calls between frozen and fixed material from each of the cases. The numbers are the number of sites average over replicates that gave concordant genotypes out of the 1494 sites on the HuSNP array.

[‡]Percentages are the average percent of concordant sites out of those with that type of call. See text for a description of how concordance was calculated.

Table 4. Number of Informative SNPs and LOH by Chromosome, Identified in the Fixed and Frozen Samples from Two Cases of Breast Cancer

Chromosome	Case 1			Case 2	
	Frozen tissue LOH [†] /informative* (DI = 1.43)	Fixed tissue		Frozen tissue LOH/informative (DI = 1.79)	Fixed tissue LOH/informative (DI = 1.76)
		LOH/informative (DI = 1.49)	LOH/informative (DI = 1.82)		
1	1/27	0/26	2/26	0/32	0/28
2	8/26	10/25	9/25	1/29	1/23
3	12/26	11/28	13/28	0/36	0/25
4	11/16	3/16	12/16	0/22	0/15
5	1/17	1/17	1/17	0/14	0/10
6	20/33	19/35	21/35	0/27	1/26
7	1/16	1/16	2/16	0/21	0/19
8	0/33	0/33	1/33	18/36	15/27
9	5/22	5/22	4/22	0/17	0/15
10	7/10	8/12	9/12	0/20	1/15
11	17/20	0/21	6/21	5/23	4/22
12	2/18	1/18	3/18	0/6	0/6
13	0/14	0/11	1/11	0/8	0/5
14	0/9	0/9	0/9	1/13	1/10
15	0/16	0/17	1/17	0/16	0/14
16	0/11	0/12	1/12	6/11	6/8
17	19/20	18/21	21/21	6/15	6/13
18	2/13	3/13	4/13	7/10	4/7
19	0/16	0/17	1/17	7/17	4/11
20	0/5	0/6	0/6	0/9	0/6
21	0/6	0/6	0/6	0/8	0/5
22	0/6	0/8	1/8	5/5	5/5
X	0/5	0/6	5/6	0/7	0/7
Unmapped	0/2	1/3	1/3	0/4	0/2
Total	106/387 (27%)	81/398 (20%)	119/398 (30%)	56/406 (14%)	45/324 (14%)

*SNP sites where the normal sample genotype call was AB (see Materials and Methods text).

[†]SNP sites which were informative in the normal and had either AA or BB genotype calls in the tumor sample (see Materials and Methods text).

with adjacent HuSNP markers, it is difficult to determine whether the apparent discordance is because of technical limitations or if the STR marker is recognizing a small region with a different allelic loss pattern than the adjacent regions scored by SNP. At four of the nine sites (D11S1394 and D17S1288 in case 1, and D17S1294 in case 2), there was at least a 5-Mb distance between the STR and SNP markers, which may be the reason for the discrepancy.

Discussion

This study examined the feasibility of using array technology, specifically the commercially available Affymetrix HuSNP array, for genome-wide allelic loss analysis of both fixed and frozen breast pathology specimens. Archival pathology specimens are a valuable resource for the genetic analysis of tumors. However, the limited quantity and quality of DNA available is a serious limitation for genetic analysis of such specimens. The quality of DNA obtained from pathology specimens is compromised by routine preservation methods that were neither designed for, nor are optimal for, DNA preservation. Formalin, the most commonly used fixative for pathology tissue specimens, has been shown to reduce the size of PCR segments that may be amplified from a sample.¹⁵ In our experience as well as in reports from the literature, DNA extracted from paraffin-embedded tissues most reliably yields PCR results in small amplicons, often under 200 nucleotides.^{15,16} Tissues frozen in OCT media for

frozen section diagnosis suffer less of a direct insult to DNA quality but are still subject to handling and storage exposures that may result in DNA fragmentation. Ideally, a genomic analysis technique for pathology specimens would maximize the data obtained from nanogram quantities of low-molecular weight DNA. Our study sought to validate array technology such as that used in the HuSNP array for use with such specimens.

In this study, samples were analyzed in duplicate to generate reliability statistics for each type of sample, and genotype data were compared between fixed and frozen samples to examine the data concordance between sample types. The HuSNP array yielded genotype results that were reliable and concordant for both fixed and frozen tumor and normal breast pathology specimens. Importantly, the DNA fragmentation that occurs with formalin fixation does not seem to affect HuSNP results, presumably because the assay relies on PCR amplicons that are shorter than 100 nucleotides in length.

In addition to analysis of genomic DNA extracted from these specimens, we also examined the data obtained from whole genome amplified material (PEP) generated from our specimens, and found similar reliability for either genomic DNA and PEP genotypes when analyzed by HuSNP. The concordance was similarly high for both genomic and PEP DNA, although slightly lower for the PEP material. This lower concordance was primarily because of an increase in no-signal genotype calls seen in the PEP material versus the genomic DNA and an indica-

tor of the potential data not obtained with amplified DNA. However, in cases in which sample is limited, using the HuSNP assay on PEP material may be an acceptable approach to genome-wide analysis. In cases in which the original sample is extremely limited, the use of whole genome amplification may make analysis possible.

We also used the data generated from the two cases to examine allelic loss in the cell populations isolated from the tumors by bivariate flow cytometry. There was an average of 379 informative SNP sites throughout the genome from all of the gDNA HuSNP assays. This is very similar to the expected distribution of heterozygosity as defined using biallelic SNP markers^{17,18} and in previously reported HuSNP data.⁴ As has been previously reported,⁴ allelic loss data obtained from the HuSNP agreed well with data obtained by the more standard method of microsatellite (STR) analysis.

Although the fixed and frozen samples from case 2 yielded highly concordant HuSNP results on all chromosomes, the cell populations with close DI (1.43 and 1.49, respectively) identified in case 1 exhibited substantial differences in LOH on chromosomes 4 and 11. Because the differences were confined to these two chromosomes, the data were not likely to be the result of a general cross-contamination, but rather reflected a biological difference between these cell populations. The second population with a DI of 1.82, identified in the formalin-fixed tissue block was also distinct and exhibited a slightly higher frequency of allelic loss throughout the genome. A diversity of cell populations is common within advanced breast tumors^{19,20} and may reflect the development of distinct genotypic clones with different behavior potential. Flow cytometric analysis can initially define cell populations with DNA content differences that can be further resolved by genomic analysis, yielding important information about tumor composition that would otherwise be obscure.

The gDNA HuSNP analysis of the two cases included in this study yielded more genome-wide data than could be obtained with a similar amount of DNA by other means, such as microsatellite marker analysis. However, it is still a low-density map, with an average of one SNP site per 8.5 Mb in the genome. Another limitation of the current HuSNP array is that many SNP sites included in the assay are clustered, so that many regions of the genome are well represented whereas others are under-represented. Given that array assays for genome-wide analyses are continuing to be developed, we expect that the next generation of genetic marker arrays using similar technology as the HuSNP will provide more uniform and higher density coverage of the genome. The data from this study indicate that future array technologies will be suitable for use with DNA obtained from routinely processed pathology specimens.

Acknowledgments

We thank Stephanie Stafford for assistance with data management and analysis, Carlo Maley for assistance in genetic map correlation, Patty Galipeau and Michael Barrett for help in the planning of the experiments included in this

study, Cassandra Neal and Jeff Delrow of the Fred Hutchinson Cancer Research Center DNA Array Facility for assistance in the hybridization and analysis of the HuSNP arrays, and Cintia De Barros for assistance with the flow sorting.

References

- Deng G, Lu Y, Zlotnikov G, Thor AD, Smith HS: Loss of heterozygosity in normal tissue adjacent to breast carcinomas. *Science* 1996, 274: 2059-2067
- Devilee P, Cornelisse CJ: Somatic genetic changes in human breast cancer. *Biochim Biophys Acta* 1994, 1198:113-130
- Osborne R, Hamshire M: A genome-wide map showing common regions of loss of heterozygosity/allelic imbalance in breast cancer. *Cancer Res* 2000, 60:3706-3712
- Lindblad-Toh K, Tanenbaum D, Daly M, Winchester E, Lui W, Villapakkam A, Stanton S, Larsson C, Hudson T, Johnson B, Lander E, Meyerson M: Loss-of-heterozygosity analysis of small-cell lung carcinomas using single-nucleotide polymorphism arrays. *Nat Biotechnol* 2000, 18:1001-1005
- Mei R, Galipeau P, Prass C, Berno A, Ghandour G, Patil N, Wolff R, Chee M, Reid B, Lockhart DJ: Genome-wide detection of allelic imbalance using human SNPs and high-density DNA arrays. *Genome Res* 2000, 10:1126-1137
- Glogovac J, Porter P, Banker D, Rabinovich P: Cytokeratin labeling of breast cancer cells extracted from paraffin embedded tissue for bivariate flow cytometric analysis. *Cytometry* 1996, 24:260-267
- Barrett MT, Reid BJ, Joslyn G: Genotypic analysis of multiple loci in somatic cells by whole genome amplification. *Nucleic Acids Res* 1995, 23:3488-3492
- Zhang L, Cui X, Schmitt K, Hubert R, Navidi W, Arnheim N: Whole genome amplification from a single cell: implications for genetic analysis. *Proc Natl Acad Sci USA* 1992, 89:5847-5851
- Porter PL, El-Bastawissi AY, Mandelson MT, Lin MG, Khalid N, Watney EA, Cousens L, White D, Taplin S, White E: Breast tumor characteristics as predictors of mammographic detection: comparison of interval- and screen-detected cancers. *J Natl Cancer Inst* 1999, 91:2020-2028
- Bauer KD, Jacobberger JW: Analysis of intracellular proteins. *Methods Cell Biol* 1994, 41:351-376
- Frank T, Svoboda-Newman S, Hsi E: Comparison of methods for extracting DNA from formalin-fixed paraffin sections for nonisotopic PCR. *Diagn Mol Pathol* 1996, 5:220-224
- Paulson T, Galipeau P, Reid B: Loss of heterozygosity analysis using whole genome amplification, cell sorting and fluorescence-based PCR. *Genome Res* 1999, 9:482-491
- Galipeau P, Cowan D, Sanchez C, Barrett M, Emond M, Levine D, Rabinovitch P, Reid B: 17p (p53) allelic losses, 4N (G2/tetraploid) populations, and progression to aneuploidy in Barrett's esophagus. *Proc Natl Acad Sci USA* 1996, 93:7081-7084
- Hahn M, Matzen S, Serth J, Pingoud A: Semiautomated quantitative detection of loss of heterozygosity in the tumor suppressor gene p53. *BioTechniques* 1995, 18:1040-1047
- Greer C, Lund J, Manos M: PCR amplification from paraffin-embedded tissues: recommendations on fixatives for long-term storage and prospective studies. *PCR Methods and Applications* 1991, 1:46-50
- Ben-Ezra J, Johnson D, Rossi J, Cook N, Wu A: Effect of fixation on the amplification of nucleic acids from paraffin-embedded material by the polymerase chain reaction. *J Histochem Cytochem* 1991, 39:351-354
- Kwok P-Y, Deng Q, Zakeri H, Taylor S, Nickerson D: Increasing the information content of STS-based genome maps: identifying polymorphisms in mapped STSs. *Genomics* 1996, 31:123-126
- Kruglyak L: The use of a genetic map of biallelic markers in linkage studies. *Nat Genet* 1997, 17:21-24
- Bonsing B, Corver W, Fleuren G, Cleton-Jansen, A-M, Devilee P, Cornelisse C: Allelotype analysis of flow-sorted breast cancer cells demonstrates genetically related diploid and aneuploid subpopulations in primary tumors and lymph node metastasis. *Genes Chromosomes Cancer* 2000, 28:173-183
- Fujii H, Marsh C, Cairns P, Sidransky D, Gabrielson E: Genetic divergence in the clonal evolution of breast cancer. *Cancer Res* 1996, 56:1493-1497